


SAFETOUCH™ IV CATHETER


**STERILE, SINGLE-USE CATHETER WITH
SAFETY DEVICE FOR IV INFUSION**

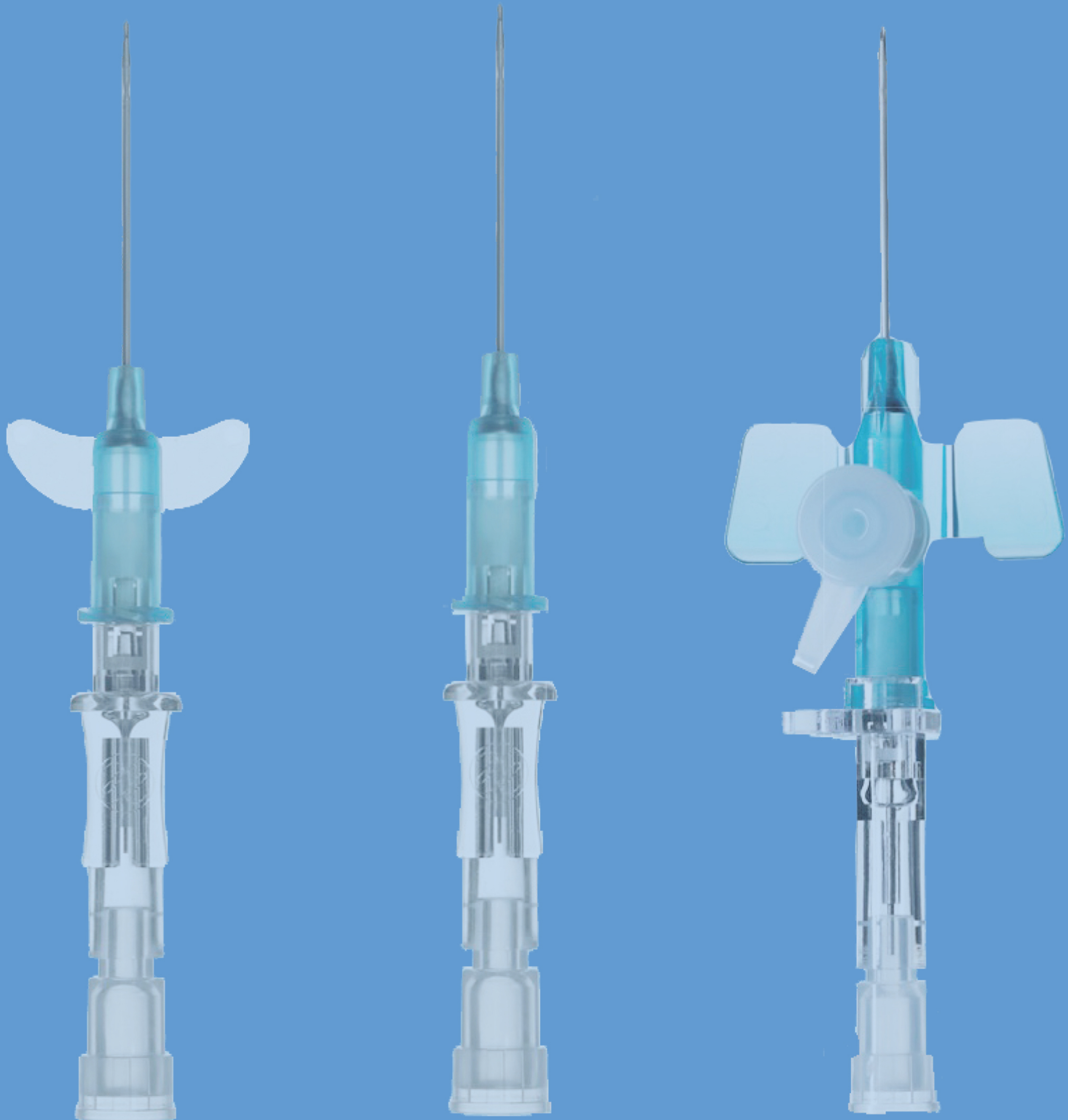

SAFETOUCH™ IV CATHETER

STERILE, SINGLE-USE CATHETER WITH SAFETY DEVICE FOR IV INFUSION

- Ultra sharp three-beveled needle, silicon coating help decrease penetration force, and make a smoothy insertion of catheter
- The canulla cannot be rotated due to its tight attachment to the catheter by a hook
- The top of the hub has a ridge that helps identify the needle facet facing up
- Transparent PUR catheter tube for quick visualization of blood in the flashback chamber to confirm the needle is in the vein
- Not made with latex, DEHP or PVC
- Irreversible, passive safety mechanism
- Audible click confirms locking of safety device
- **More secure with two-layer protectors (metal layer on the inside, plastic player on the outside)**
- **Limit cross-contamination thanks to the plastic protective layer that prevents blood from spilling**
- 3-stripe radiopaque, which is visualized via X-ray image in case the fractured tips are left in the body
- Available in 16-24G, The hub color helps to identify the size of the catheter
- IV infusion duration: hours to several days
- EO sterile, individually packaged
- CE, ISO 13485 certificates
- Due date: 5 years

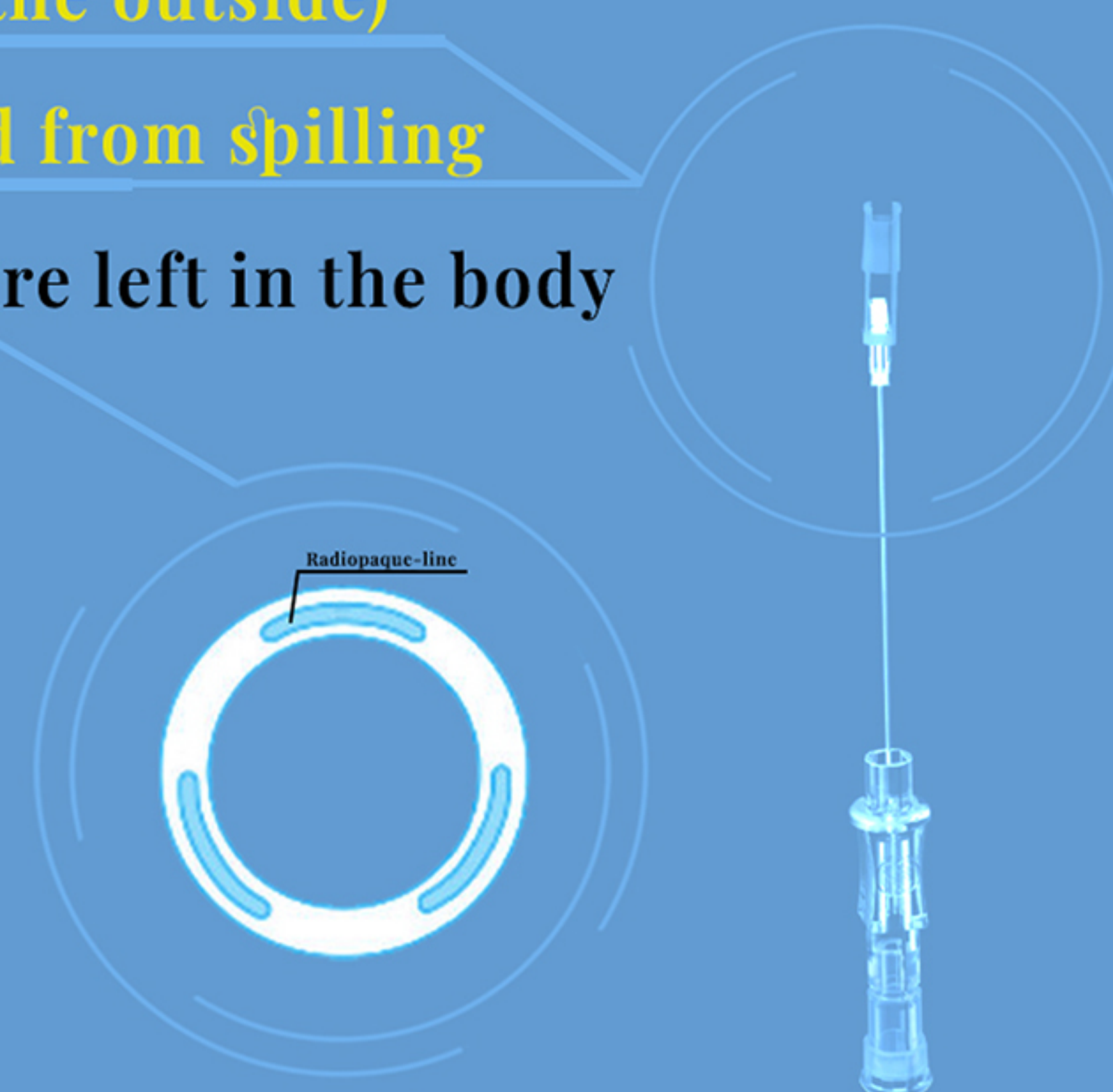

PRODUCT RANGE OVERVIEW

Code	Description	Flow rate not less than	Color
Safetouch Cath (Straight)			
C-SNIC-1650-ECE	SAFETOUCH CATH (STRAIGHT) 16Gx2" (50mm), ETO (CE)	≥ 162 ml/phút	
C-SNIC-1832-ECE	SAFETOUCH CATH (STRAIGHT) 18Gx1-1/4" (32mm), ETO (CE)	≥ 81 ml/phút	
C-SNIC-2032-ECE	SAFETOUCH CATH (STRAIGHT) 20Gx1-1/4" (32mm), ETO (CE)	≥ 49 ml/phút	
C-SNIC-2225-ECE	SAFETOUCH CATH (STRAIGHT) 22Gx1" (25mm), ETO (CE)	≥ 29 ml/phút	
C-SNIC-2232-ECE	SAFETOUCH CATH (STRAIGHT) 22Gx1-1/4" (32mm), ETO (CE)	≥ 29 ml/phút	
C-SNIC-2419-ECE	SAFETOUCH CATH (STRAIGHT) 24Gx3/4" (19mm), ETO (CE)	≥ 17 ml/phút	
C-SNIC-2425-ECE	SAFETOUCH CATH (STRAIGHT) 24Gx1" (25mm), ETO (CE)	≥ 17 ml/phút	
Safetouch Cath wing cath without injection port			
C-SWC-1650-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 16Gx2" (50mm), ETO (CE)	≥ 162 ml/phút	
C-SWC-1832-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 18Gx1-1/4" (32mm), ETO (CE)	≥ 81 ml/phút	
C-SWC-1850-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 18Gx2" (50mm), ETO (CE)	≥ 81 ml/phút	
C-SWC-2032-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 20Gx1-1/4" (32mm), ETO (CE)	≥ 49 ml/phút	
C-SWC-2050-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 20Gx2" (50mm), ETO (CE)	≥ 49 ml/phút	
C-SWC-2225-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 22Gx1" (25mm), ETO (CE)	≥ 29 ml/phút	
C-SWC-2232-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 22Gx1-1/4" (32mm), ETO (CE)	≥ 29 ml/phút	
C-SWC-2419-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 24Gx3/4" (19mm), ETO (CE)	≥ 17 ml/phút	
C-SWC-2425-ECE	SAFETOUCH WING CATH W/O INJECTION PORT 24Gx1" (25mm), ETO (CE)	≥ 17 ml/phút	
Safetouch Cath wing with injection port			
C-SWCP-1650-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 16Gx2" (50mm), ETO (CE)	≥ 162 ml/phút	
C-SWCP-1832-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 18Gx1-1/4" (32mm), ETO (CE)	≥ 81 ml/phút	
C-SWCP-1850-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 18Gx2" (50mm), ETO (CE)	≥ 81 ml/phút	
C-SWCP-2032-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 20Gx1-1/4" (32mm), ETO (CE)	≥ 49 ml/phút	
C-SWCP-2050-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 20Gx2" (50mm), ETO (CE)	≥ 49 ml/phút	
C-SWCP-2225-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 22Gx1" (25mm), ETO (CE)	≥ 29 ml/phút	
C-SWCP-2232-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 22Gx1-1/4" (32mm), ETO (CE)	≥ 29 ml/phút	
C-SWCP-2419-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 24Gx3/4" (19mm), ETO (CE)	≥ 17 ml/phút	
C-SWCP-2425-ECE	SAFETOUCH WING CATH WITH INJECTION PORT 24Gx1" (25mm), ETO (CE)	≥ 17 ml/phút	

MATERIALS USED

Needle (cannula)

Material: Stainless steel SUS-304

Cutting angle

Three-beveled

Lubricant

Material: Silicone Polydimethylsiloxane

Catheter tube

Material: Polyurethane (PUR)

Caulking pin

Material: Stainless steel SUS-304

Ledge

Material: Polypropylene (PP)

Safety ring

Material: Polyoxymethylene (POM)

Catheter hub

Material: Polypropylene (PP)

Hook

Material: Polycarbonate (PC)

Protector

Material: Stainless steel SUS-304

Needle hub


Material: Polycarbonate (PC)

**Microglass filter
prevents blood to leak out**


Material: Random microglass filter

Vent fitting

Material: Polypropylene (PP)

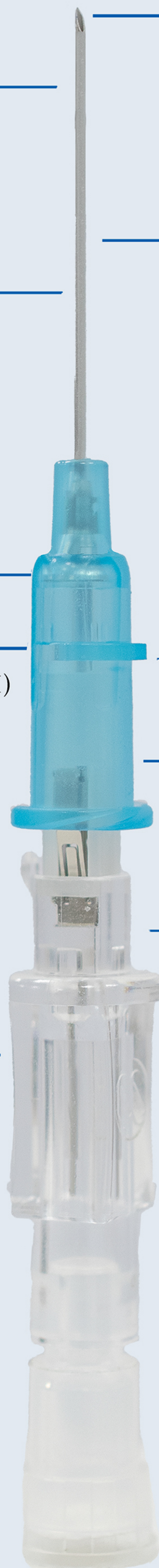

**AmSafe - Pre-filled
Normal Saline Flush Syringe
0.9% Sodium Chloride Injection**

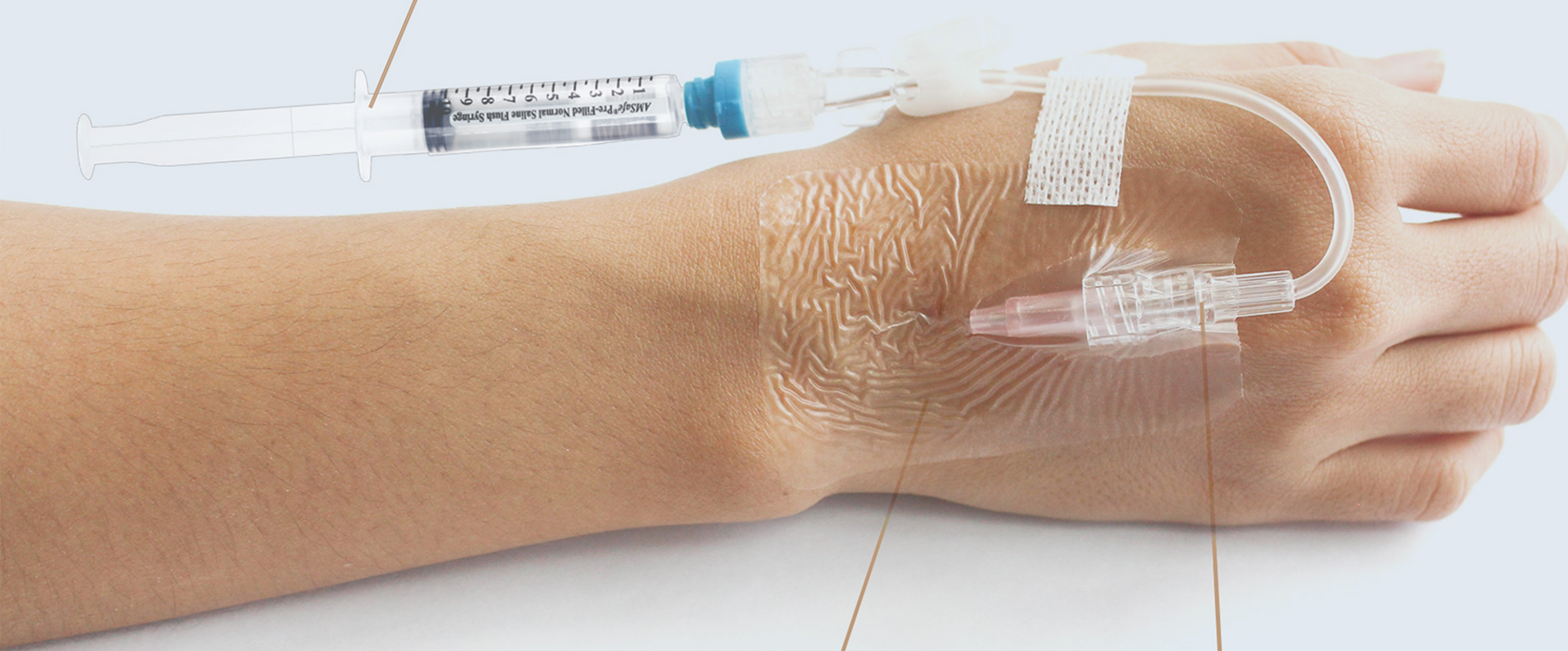

Safe derm

**Safetouch plug
with extension tube**



NIPRO CORPORATION

3-9-3 Honjo-nishi, Kita-ku, Osaka 531-8510, Japan
Tel: +81-6-6375-6713

NIPRO ASIA LTD., PTE

8 Temasek Boulevard #40-02 Suntec Tower 3
Singapore 038988
Tel: +65 6588-2281

NIPRO SALES VIETNAM CO., LTD

Floor 15th, Viettel Complex Bulding,
285 Cach Mang Thang Tam, Ward 12, District 10, HCMC
Tel: +028 6277 0057